

io

March 14, 2017

```
In [ ]: %gui qt
```

1 Required knowledges

- numpy (basic knowledges)
- Qt (basic knowledges)
- h5py (basic knowlegdes)

2 Silx IO API

```
In [ ]: import silx.io
```

2.1 Open a file

```
In [ ]: # open
obj = silx.io.open("data/test.h5")
# do your stuff here
obj.close()
```

2.2 Open a file with context manager

```
In [ ]: # or using context manager
with silx.io.open("data/test.h5") as obj:
    # do your stuff here
    # the close is called for you at the end of the with
pass
```

2.3 Common properties

```
In [ ]: obj = silx.io.open("data/test.h5")

obj.name    # path from file to the object
obj.parent  # hdf5 group containing the object
obj.file    # hdf5 file containing the object

# test object type
if silx.io.is_file(obj):
```

```

    # this is a root file
    # path of the file from the file system
    obj.filename

    if silx.io.is_group(obj):
        # this is a group
        # BTW a file is a group
        pass

    if silx.io.is_dataset(obj):
        # this is a dataset\
        pass

```

2.4 Node traversal

```

In [ ]: obj = silx.io.open("data/test.h5")

    if silx.io.is_group(obj):
        # it contains child

        # number of child
        len(obj)

        # iterator on child names
        obj.keys()

        # access to a child
        child = obj["arrays"]

        # access to a child using a path
        child = obj["arrays/float_3d"]

        # the path can be absolute
        child = obj["/arrays/float_3d"]

```

2.5 Data access

```

In [ ]: h5 = silx.io.open("data/test.h5")
        obj = h5["arrays/float_3d"]

    if silx.io.is_dataset(obj):
        # it contains data

        # a dataset provides information to the data
        obj.shape    # multidimensional shape
        obj.size     # amount of items
        obj.dtype    # type of the array

```

```

    # copy the full data as numpy array
    data = obj[...]

    # or a part of it (using numpy selector)
    data = obj[1:2, ::3, 2]

    scalar = h5["scalars/int64"]
    if silx.io.is_dataset(scalar):
        # scalar dataset have an empty shape
        assert scalar.shape == ()

    # special case to access to the value of a scalar
    data = scalar[()]

```

2.6 Spec file as HDF5

```

In [ ]: import silx.io

h5like = silx.io.open('data/oleg.dat')

# print available scans
print(h5like['/'].keys())

# print available measurements from the scan 94.1
print(h5like['/94.1/measurement'].keys())

# get data from measurement
time = h5like['/94.1/measurement/Epoch']
bpm = h5like['/94.1/measurement/bpmi']
mca = h5like['/94.1/measurement/mca_0/data']

```

2.7 EDF file as HDF5

```

In [ ]: import silx.io.utils

h5like = silx.io.open("data/ID16B_diatomee.edf")

# here is the data as a cube using numpy array
# it's a cube of images * number of frames
data = h5like["/scan_0/instrument/detector_0/data"]
# here is the first image
data[0]
print(data[0].shape)

# groups containing datasets of motors, counters
# and others metadata from the EDF header
motors = h5like["/scan_0/instrument/positioners"]
counters = h5like["/scan_0/measurement"]

```


```

others    = h5like[("/scan_0/instrument/detector_0/others")]
print("motor names", list(motors.keys()))

# reach a motor named 'srot'
# it's a vector of values * number of frames
srot = motors["srot"]
# here is the monitor value at the first frame
srot[0]
print(srot[...])

```

3 Silx HDF5 tree

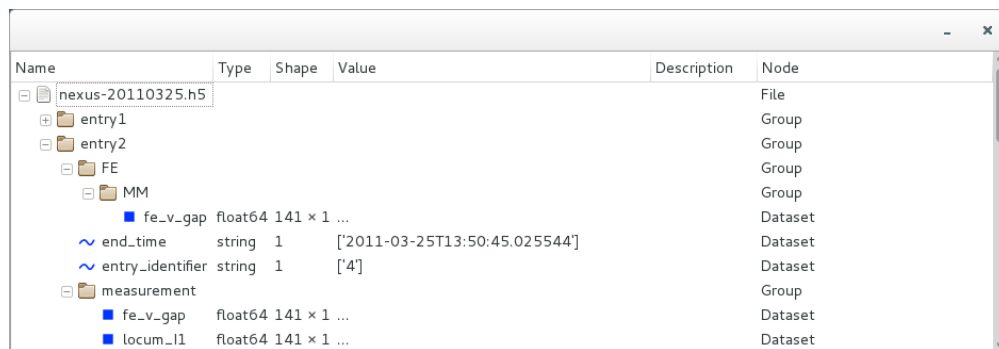

Name	Type	Shape	Value	Description	Node
nexus-20110325.h5					File
entry1					Group
entry2					Group
FE					Group
MM					Group
fe_v_gap	float64	141 × 1 ...			Dataset
end_time	string	1	['2011-03-25T13:50:45.025544']		Dataset
entry_identifier	string	1	['4']		Dataset
measurement					Group
fe_v_gap	float64	141 × 1 ...			Dataset
locum_l1	float64	141 × 1 ...			Dataset

HDF5 tree screenshot

- Getting start with the HDF5 tree (http://pythonhosted.org/silx/modules/gui/hdf5/getting_started.html)

3.1 Create the widget

```

In [ ]: import silx.gui.hdf5
        tree = silx.gui.hdf5.Hdf5TreeView()
        tree.setVisible(True)
        model = tree.findHdf5TreeModel()

```

3.2 Feed it with an HDF5

```

In [ ]: h5 = silx.io.open("data/test.h5")
        model.insertH5pyObject(h5)

```

3.3 Feed it with a Spec file

```


In [ ]: h5 = silx.io.open("data/oleg.dat")
        model.insertH5pyObject(h5)

```

3.4 Feed it with an EDF file

```
In [ ]: h5 = silx.io.open("data/ID16B_diatomee.edf")
        model.insertH5pyObject(h5)
```

4 Silx DataViewer

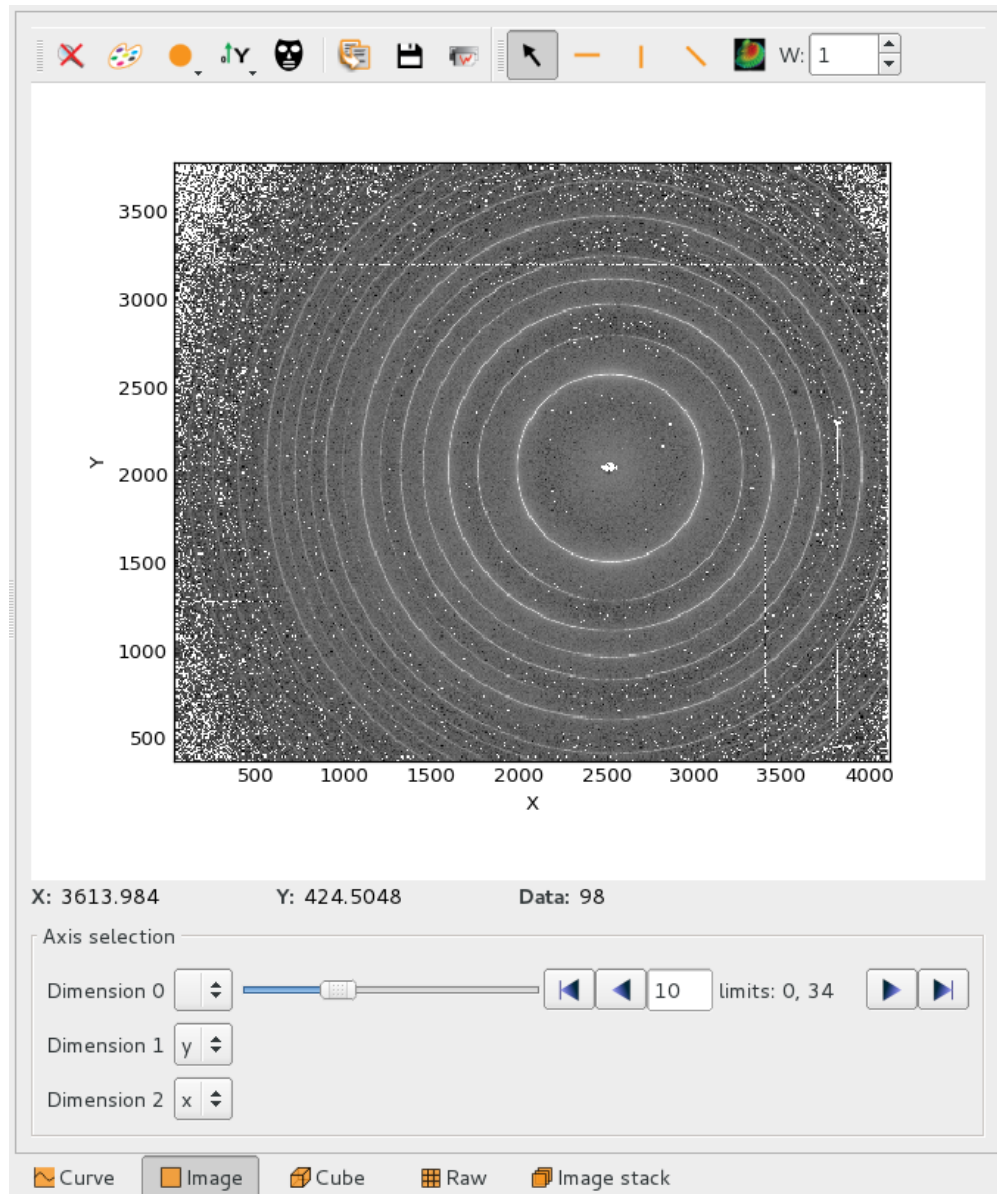

DataViewer screenshot

4.1 Create the widget

```
In [ ]: import silx.gui.data.DataViewerFrame
        dataviewer = silx.gui.data.DataViewerFrame.DataViewerFrame()
        dataviewer.setVisible(True)
        dataviewer.resize(500, 500)
```

4.2 Feed it with a numpy array

```
In [ ]: import numpy
        data = numpy.random.rand(100, 100, 100)
        dataviewer.setData(data)
```

4.3 Feed it with a HDF5 dataset

```
In [ ]: import silx.io
        h5like = silx.io.open("data/ID16B_diatomee.h5")
        dataset = h5like["/data/0299"]
        dataviewer.setData(dataset)
```

5 Exercises

The exercise is based on a phase contrast data. It will help you to create a custom application to browse data.

5.1 Exercise 1

- Browse an HDF5 file
- Use [getting started with HDF5 widgets](#)
- Identify path of the data
- Access to the data

```
In [ ]: import silx.io

        #
        # EXERCISE: Open the file 'data/ID16B_diatomee.h5'
        #

        h5 = ...

        #
        # EXERCISE: Display the file into the HDF5 tree
        #

        from silx.gui import hdf5
        tree = hdf5.Hdf5TreeView()
        model = tree.findHdf5TreeModel()
        ...
        tree.setVisible(True)
```

```

In [ ]: #
        # EXERCISE: Access to one frame of the image
        #

        print(...)

        #
        # EXERCISE: Display it with the data viewer
        #
        import silx.gui.data.DataViewerFrame
        viewer = silx.gui.data.DataViewerFrame.DataViewerFrame()
        ...
        viewer.setVisible(True)

```

5.2 Exercise 2

1. From the HDF5 tree, identify path name for
 - one data frame
 - one background
 - one flatfield
2. Compute flatfield correction

The computation of corrected images is done using this equation using data, flatfield, and background information.

$$corrected = \frac{data - background}{flatfield - background}$$

```

In [ ]: def correctedImage(data, background, flatfield):
        data = numpy.array(data, dtype=numpy.float32)
        flatfield = numpy.array(flatfield, dtype=numpy.float32)
        return (data - background) / (flatfield - background)

        #
        # EXERCISE: Reach one data frame, a background and a flatfield from 'data/'
        #

        ...

        #
        # EXERCISE: Compute the corrected image
        #

        ...

        #
        # EXERCISE: Display it with the data viewer
        #

```

...

5.3 Exercise 3

1. Connect together an HDF5 tree view and a data viewer
 - Use [getting started with HDF5 widgets](#)

```
In [ ]: from silx.gui import qt
        from silx.gui import hdf5
        import silx.gui.data.DataViewerFrame

class ViewerEx3(qt.QMainWindow):

    def __init__(self, parent=None):
        qt.QMainWindow.__init__(self, parent)
        widget = self.createCentralWidget()
        self.setCentralWidget(widget)

    def createCentralWidget(self):
        splitter = qt.QSplitter(self)

        # the tree
        self.tree = silx.gui.hdf5.Hdf5TreeView(self)
        # the data viewer
        self.viewer = silx.gui.data.DataViewerFrame.DataViewerFrame(self)

        splitter.addWidget(self.tree)
        splitter.addWidget(self.viewer)
        splitter.setStretchFactor(1, 1)

        #
        # EXERCISE: Connect the callback onTreeActivated (bellow)
        #               to a mouse event from the tree
        #

        return splitter

    def onTreeActivated(self):

        #
        # EXERCISE: Reach selected objects from the tree
        #

        #
        # EXERCISE: Provide it to the data viewer
        #
```



```

        pass

    def appendFile(self, filename):
        model = self.tree.findHdf5TreeModel()
        model.insertFile(filename)
        print("Load %s" % filename)

In [ ]: viewer = ViewerEx3()
        viewer.appendFile('data/ID16B_diatomee.h5')
        viewer.setVisible(True)

```

5.4 Exercise 4

1. Use the previous application to display corrected data

```

In [ ]: class ViewerEx4(ViewerEx3):

    def onTreeActivated(self):
        selectedObjects = list(self.tree.selectedH5Nodes())
        if len(selectedObjects) == 0:
            self.viewer.setData("Nothing selected")

        elif len(selectedObjects) > 1:
            self.viewer.setData("Too much things selected")

        else:
            obj = selectedObjects[0]
            node = obj.h5py_object

            if "/data/" in node.name:
                # That's a data from the /data group
                data = self.computeCorrectedImage(node)
                self.viewer.setData(data)
            else:
                # Other data is displayed in a normal way
                self.viewer.setData(obj)

    def computeCorrectedImage(self, h5data):
        """
        :param h5data: H5py dataset selected from the group /data/
        """
        background = self.getBackground(h5data)
        flatfield = self.getFlatField(h5data)

        raw = numpy.array(h5data, dtype=numpy.float32)
        flatfield = numpy.array(flatfield, dtype=numpy.float32)
        background = background[...]
        return (raw - background) / (flatfield - background)

```

```

def getBackground(self, h5data):
    """
    :param h5data: H5py dataset selected from the group /data/
    """

    #
    # EXERCISE: Return the background image from the dataset
    #

    return None

def getFlatField(self, h5data):
    """
    :param h5data: H5py dataset selected from the group /data/
    """

    #
    # EXERCISE: Return the flatfield image from the dataset
    #
    #
    # 1) you can return a flatfield by default
    # 2) you can return the closest flat field according to t
    # 3) you can return an interpolation of the 2 flatfields

    return None

```

```

In [ ]: viewer = ViewerEx4()
viewer.appendFile('data/ID16B_diatomee.h5')
viewer.setVisible(True)

```